



A Publication
of Reliable Methods
for the Preparation
of Organic Compounds

Working with Hazardous Chemicals

The procedures in *Organic Syntheses* are intended for use only by persons with proper training in experimental organic chemistry. All hazardous materials should be handled using the standard procedures for work with chemicals described in references such as "Prudent Practices in the Laboratory" (The National Academies Press, Washington, D.C., 2011; the full text can be accessed free of charge at http://www.nap.edu/catalog.php?record_id=12654). All chemical waste should be disposed of in accordance with local regulations. For general guidelines for the management of chemical waste, see Chapter 8 of Prudent Practices.

In some articles in *Organic Syntheses*, chemical-specific hazards are highlighted in red "Caution Notes" within a procedure. It is important to recognize that the absence of a caution note does not imply that no significant hazards are associated with the chemicals involved in that procedure. Prior to performing a reaction, a thorough risk assessment should be carried out that includes a review of the potential hazards associated with each chemical and experimental operation on the scale that is planned for the procedure. Guidelines for carrying out a risk assessment and for analyzing the hazards associated with chemicals can be found in Chapter 4 of Prudent Practices.


The procedures described in *Organic Syntheses* are provided as published and are conducted at one's own risk. *Organic Syntheses, Inc.*, its Editors, and its Board of Directors do not warrant or guarantee the safety of individuals using these procedures and hereby disclaim any liability for any injuries or damages claimed to have resulted from or related in any way to the procedures herein.

These paragraphs were added in September 2014. The statements above do not supersede any specific hazard caution notes and safety instructions included in the procedure.

Organic Syntheses, Coll. Vol. 3, p.608 (1955); Vol. 27, p.56 (1947).

N-METHYL-1-NAPHTHYLCYANAMIDE

[Cyanamide, methyl-(1-naphthyl)-]

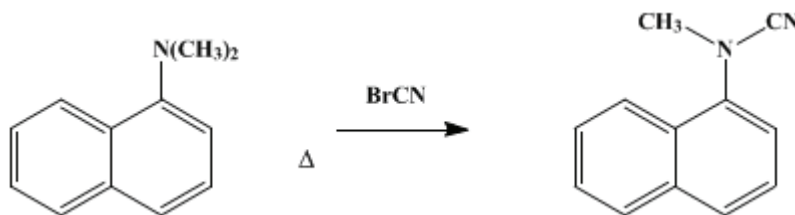

Submitted by Homer W. J. Cressman

Checked by Nathan L. Drake and Werner R. Boehme..

1. Procedure

Caution! Cyanogen bromide is highly toxic. This preparation should be carried out under a good hood.

A mixture of 171 g. (1 mole) of dimethyl- α -naphthylamine (Note 1) and 125 g. (1.2 moles) of cyanogen bromide¹ in a 1-l. flask (Note 2) is heated under reflux on the steam bath for 16 hours (Note 3). The cooled reaction mixture is added to 2.5 l. of dry ether, and the insoluble quaternary salt (Note 4) is filtered. The ether filtrate is extracted with four 800-ml. portions of approximately 15% hydrochloric acid solution (Note 5) and washed with five 500-ml. portions of water. The ether solution is dried with 30–35 g. of anhydrous calcium sulfate (Drierite). After filtration, the solvent is removed by distillation at atmospheric pressure on the steam bath; the residue is fractionated under reduced pressure. The yield of pale yellow oil, boiling at 170–171°/1 mm. (185–187°/3 mm.), amounts to 115–122 g. (63–67%) (Note 6).

2. Notes

- Better yields are obtained by dealkylating the tertiary amine in this manner than by starting with the secondary amine. In a similar manner, *N*-ethyl-1-naphthylcyanamide, b.p. 165–168°/2 mm., is obtained in a 48% yield from diethyl- α -naphthylamine; the yield from the secondary amine is only 25%.
- Ground-glass equipment is preferable.
- In order to avoid loss of cyanogen bromide an efficient condenser must be used and the heating should be very gradual at first.
- About 10 g. of α -naphthyltrimethylammonium bromide, m.p. 160°, is generally obtained.
- About 25 g. of crystalline dimethyl- α -naphthylamine hydrochloride may be recovered by efficient chilling of the hydrochloric acid extracts.
- The *N*-methyl-1-naphthylcyanamide turns dark green upon standing exposed to the air for several days. It may be preserved in sealed, evacuated ampoules in the absence of light.

3. Discussion

The preparation described is based on the method of von Braun and co-workers.²

This preparation is referenced from:

- *Org. Syn. Coll. Vol. 3, 609*

1. *Org. Syntheses Coll. Vol. 2*, 150 (1943).
 2. Von Braun, Heider, and Muller, *Ber.*, **51**, 281 (1918).
-

Appendix
Chemical Abstracts Nomenclature (Collective Index Number);
(Registry Number)

dimethyl- α -naphthylamine

diethyl- α -naphthylamine

dimethyl- α -naphthylamine hydrochloride

hydrochloric acid (7647-01-0)

ether (60-29-7)

calcium sulfate (7778-18-9)

Cyanogen bromide (506-68-3)

N-Methyl-1-naphthylcyanamide,
Cyanamide, methyl-(1-naphthyl)- (53663-33-5)

N-ethyl-1-naphthylcyanamide

α -naphthyltrimethylammonium bromide