



GC Conditions:

Chiral BETA DEX™ 325 fused silica capillary column (30 m x 0.25 mm x 0.25 mm film thickness). Temperature gradient: 3 minutes of equilibration at 25 °C, followed by linear increase to 70 °C over 90 minutes (0.5 °C/min).

Scalemic Mixture

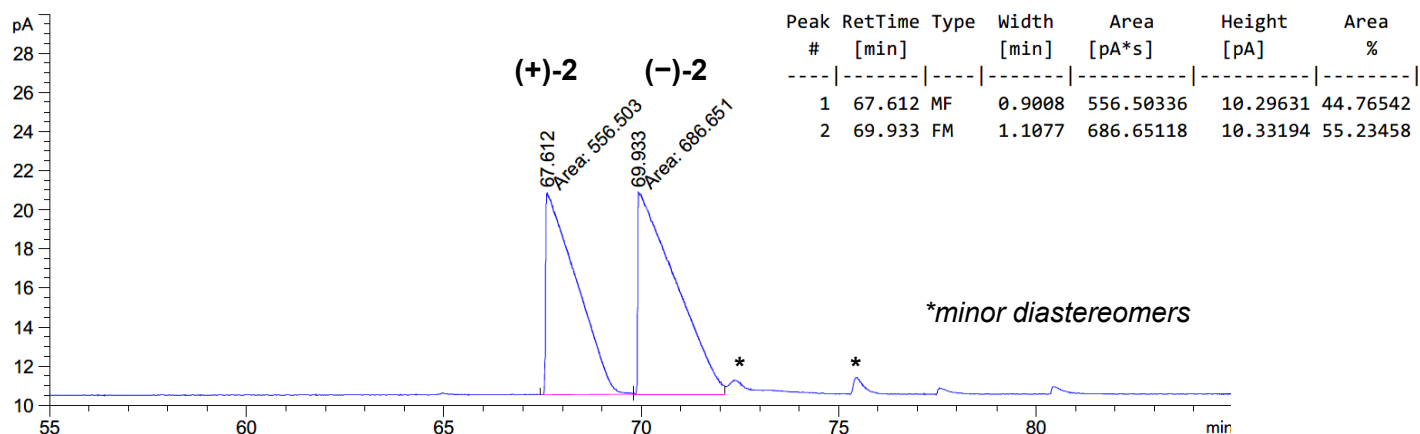

(+)-2 (S)-1-((R)-oxiran-2-yl)prop-2-en-1-ol

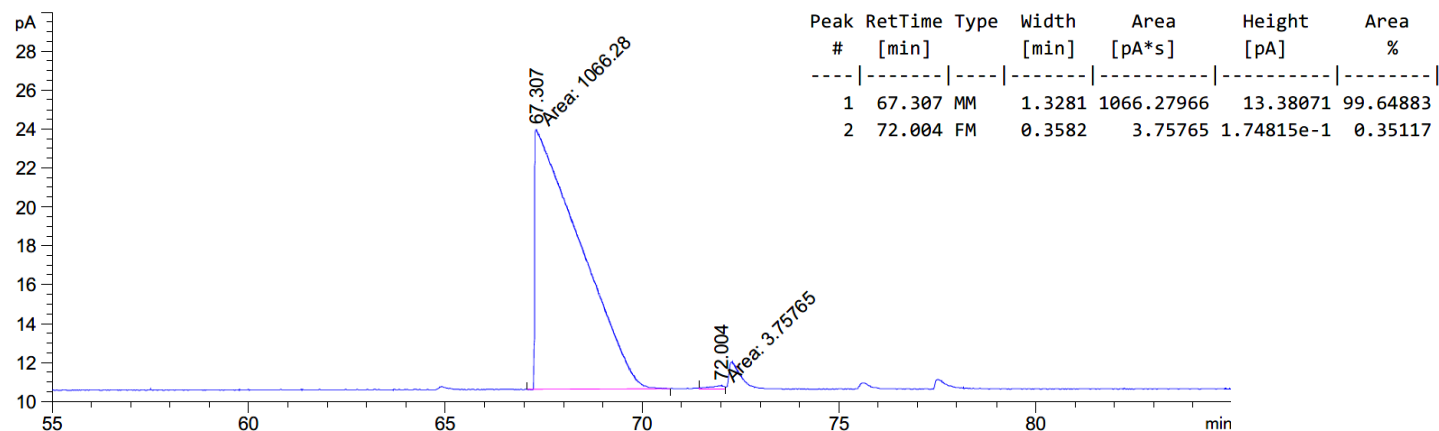

(-)-2 (R)-1-((S)-oxiran-2-yl)prop-2-en-1-ol

